

Coronis Fusion 6MP (MDCC-6530)

6MP diagnostic color display system



- **Color fusion display for flexible multimodality imaging**
- **With toolset to spark productivity**
- **Automated QA and calibration**

A brighter diagnosis

Coronis Fusion, Barco's renowned multi-modality display for radiologists, now comes in a new, energy-efficient and lightweight design. Bright on so many levels, Coronis Fusion has been designed to help radiologists provide care with confidence.

See more details quickly thanks to the high brightness, high contrast ratio and best-in-class image quality. What's more, the wide color gamut in combination with SteadyColor™ calibration technology helps you see even more colors and more details on the 30-inch screen.

Smart workflow


Designed for your comfort and to boost productivity, Coronis Fusion comes with smart image-enhancing features and workflow tools. The display's wide viewing angle combined with the SoftGlow™ task and wall light help reduce eye strain. Thanks to SpotView, radiologists can further improve detection accuracy as well as reading productivity.

Effortless quality & compliance

Like all of Barco's medical display systems, Coronis Fusion comes with QAWeb Enterprise, a cloud-based technology for automated calibration, Quality Assurance and compliance to ensure maximum uptime of the display with no need for human intervention.

A one-stop-shop solution

Barco medical display systems include exclusive display controllers that are validated with the latest workstations and with all major PACS applications. All components are under the full 5-year warranty (including the display backlight) for complete peace of mind.

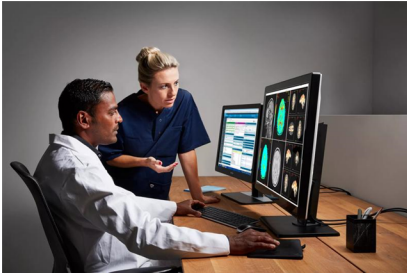

Ensuring diagnostic confidence with MDR Class IIa

Our radiology displays are MDR-certified as Class IIa. Their product information has been reviewed and cleared by independent medical and technical experts, and is audited yearly. In other words, we ensure diagnostic confidence and peace of mind for our users.

Models

- MDCC-6530: includes protective, anti-reflective glass cover
- MDCC-6530 NC: without cover

Product specifications**CORONIS FUSION 6MP (MDCC-6530)****General specifications**

Screen technology	IPS
Active screen size (diagonal)	772 mm (30.4")
Active screen size (H x V)	654 x 409 mm (25.8 x 16.1")
Aspect ratio (H:V)	16:10
Resolution	Native 6MP (3280 x 2048 pixels) Configurable to 2 x 3MP+ (1640 x 2048 pixels) Configurable to 2 x 3MP (1536 x 2048 pixels)
Pixel pitch	0.1995 mm
Color imaging	Yes
Gray imaging	Yes
Bit depth	30 bit
Viewing angle (H, V)	178°
Uniformity correction	Color PPU
SteadyColor	Yes, when used as system with MXRT-x600 display controller or later, using Barco MXRT driver 10.171.1.0 or later and with QAWeb Enterprise 2.2 or later.
Ambient light presets	Yes, reading room selection
Ambient light sensor	Yes
Front sensor	Yes, I-Guard
Maximum luminance (panel typical)	1050 cd/m ²
DICOM calibrated luminance	600 cd/m ²
Contrast ratio (panel typical)	2000:1
Response time ((Tr + Tf)/2) (typical)	18 ms
Housing color	Black / White
Video input signals	2x DisplayPort 1.2
Video output signals	1x DisplayPort (MST)
USB ports	1x USB 2.0 upstream (endpoint) 2x USB 2.0 downstream 1x USB 2.0 downstream with high-power charging functionality
Power rating	100-240 Vac, 50/60 Hz, 3.6-1.6 A
Power consumption	80 W (nominal) @ calibrated luminance of 600 cd/m ² < 0.5 W (hibernate) < 0.5 W (standby)
Dimensions with stand (W x H x D)	714 x 524~624 x 240 mm
Dimensions w/o stand (W x H x D)	714 x 478 x 74 mm
Dimensions packaged (W x H x D)	800 x 650 x 295 mm
Net weight with stand	MDCC-6530: 17.7 kg MDCC-6530 NC: 16.3 kg
Net weight w/o stand	MDCC-6530: 13.1 kg MDCC-6530 NC: 11.7 kg
Net weight packaged	MDCC-6530: 22.3 kg (without optional accessories) MDCC-6530 NC: 20.9 kg (without optional accessories)
Tilt	-5° to +25°
Swivel	-30° to +30°
Pivot	N/A
Height adjustment range	100 mm
Mounting standard	VESA (100 mm)
Screen protection	MDCC-6530: Protective, anti-reflective glass cover MDCC-6530 NC: No glass cover
Recommended modalities	All digital images, except digital mammography

Product specifications

CORONIS FUSION 6MP (MDCC-6530)

Certifications	<div>FDA 510(k) K191845 CE0123 (Medical Device) CCC (China) INMETRO (Brazil -Prodcut numbers K9301636B(1X), K9301638B(1X) KC (Korea) BIS (India) EAC (Russia, Kazakhstan, Belarus, Armenia and Kyrgyzstan) Safety specific: IEC 60950-1:2005 + A1:2009+A2:2013 EN 60950-1:2006 + A1:2010 + A11:2009 + A12:2011 + A2:2013 IEC 62368-1:2018 EN IEC 62368-1:2020 + A11:2020 IEC 60601-1:2005+A1:2012+A2:2020 EN 60601-1:2006 + A1:2013 + A12:2014+A2:2021 AAMI ES 60601-1:2005 + A1:2012+A2:2021 CAN/CSA C22.2 No. 60601-1:2014 (Reaffirmed 2022) EMI specific: IEC 60601-1-2:2014+A1:2020 (ed.4.1) EN 60601-1-2:2015+A1:2021 (ed.4.1) FCC part 15 Class B ICES-001 Level B VCCI (Japan) Environmental: EU RoHS China RoHS REACH Canada Health WEEE Packaging Directive</div>
Supplied accessories	<div>User guide Documentation disc System sheet Video cables Mains cables USB cable</div>
Optional accessories	<div>Graphics board Touch pad</div>
QA software	<div>QAWeb</div>
Warranty	<div>5 years, including 40000 hrs backlight warranty</div>
Operating temperature	<div>0 °C to 35 °C (10 °C to 30 °C within specs)</div>
Storage temperature	<div>-20 °C to 60 °C</div>
Operating humidity	<div>20% to 85% (non-condensing)</div>
Storage humidity	<div>20% to 85% (non-condensing)</div>
Operating pressure	<div>70 kPa minimum</div>
Storage pressure	<div>50 to 106 kPa</div>

Last updated: 08 Jul 2024

© 2024 Barco nv. All rights reserved. Reproduction in whole or in part without written permission is prohibited. All brand names and product names are trademarks, registered trademarks or tradenames of their respective holders. Due to continued innovation, information and technical specifications are subject to change without prior notice. Please check www.barco.com for the latest specifications.

